

Carrot 3000

SAFETY FEATURES

Our most comprehensive special needs car seat.

Ergonomically designed **HEAD SUPPORT** holds the head and protects from side impacts and centrifugal forces

SHOULDER PROTECTOR WINGS: Extendable wings hold and protect the body from side impacts and centrifugal forces

TETHER OPTIONS include additional security tether for the top of backrest or ISOfix latch connectors. These may be required by your local authorities

Large **CHEST SUPPORT** holds the body securely

2 positions of **SEAT DEPTH ADJUSTMENT** for a better fit

RECLINE FUNCTION with angle adjustable seat and backrest



Carrot 3000

The Carrot 3000 caters for safety and the importance of keeping the child in a supported posture for travel.

- Ensure Carrot 3000 **back is in contact with the vehicle seat back** (no big gap behind headrest).
- Strap at the rear of the Carrot 3000's **headrest is to be secured around vehicle's headrest**, or backrest.
- If the car has individual separation in the rear seats, the **lower horizontal strap at the rear of the Carrot 3000 is to be secured**.
- Ensure the **shoulder straps are firmly secure**, and are level with the child's shoulder height, if these are below the child's shoulder height, the seat needs to be grown – Please contact your Occupational Therapy team.
- If a cross strap is being used, please ensure this is **sitting no higher than the top being level with the child's underarm**.
- Ensure child has **not outgrown the weight capacity which is 15-36kg**. If child is approaching this weight limit, consider what car seat would be next option.
- Ensure the **car's seatbelt is used EVERY time**. This is to be **securely fastened OVER the child and seat**. Please note this seat does not have capacity to have a belt threaded through the back as it is a booster seat. You are able to purchase and install a padded sleeve for the seatbelt if you feel the belt is going to irritate your child.
- Please **note the red belt guides are to be used to ensure the belt is sitting correctly**, is kept low across the child's hips and is kept off the child's face.
- If this unit is fitted with a **swivel base check the latch is working correctly**, and that the car seat swivels freely.
- Periodically **check for loose or missing parts**.


Contact the car seat specialists for more information:
0800 543 343
solutions@medifab.com

Distributed by **medifab**[®]

The quick reference guide does not replace the requirement to read the user manual and safety information therein.



Ensure back is in contact with the vehicle seat back

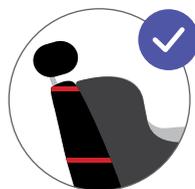

Ensure back straps are firmly fixed around the vehicle's headrest and seat



Ensure shoulder straps are firmly secure



Ensure child has not outgrown the weight capacity limit



The car's seatbelt must be worn over top of the child and car seat harness